



The logo for Medicanix, featuring the word "MEDICANIX" in a bold, white, sans-serif font. A red and orange curved line arches over the text, starting under the 'M' and ending under the 'X'.
MEDICANIX

Your Veterinary Equipment Specialist




Scope of our Services:

Medicanix has offered New York, New Jersey and Connecticut with Veterinary equipment repair and service since 1982. Many of our clients have high equipment demands which Medicanix strives to meet the needs of.

Veterinarians use a large variety of equipment and much of it gets heavy usage and excessive wear. Veterinarians can experience problems with equipment such as anesthesia systems / vaporizers that periodically need calibration, stainless veterinary tables at times leak fluid and/or drift down, veterinary ceiling mounted surgery lights burn out and break down. The autoclave in a busy veterinary practice often gets a lot of usage and requires periodic servicing and the lab microscope can often get overly dirty needing a good cleaning. The dog and cat scales need to be periodically calibrated, vital signs monitors often need new cables/lead-wires, oxygen hoses / fittings break and start leaking and the list goes on and never ends. We understand that you need to run a successful practice and let us help you with this ongoing list of equipment nightmares so they do not drag your practice down.

Medicanix repairs it all when it comes to veterinary equipment. Our comprehensive loaner fleet can keep your practice up and running during any equipment malfunction.

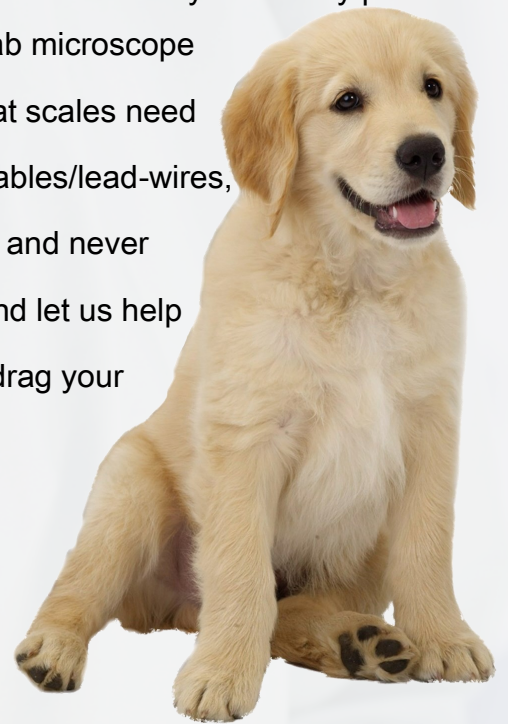

OUR GUARANTEE:

Is our commitment to our customers, that they are 100% satisfied and the workmanship rendered is second to none.

Three Equipment Maintenance Services for \$ 399

One Service Visit - One Price (price is for labor service only - replacement parts are additional)



Anesthesia System
Calibration Checks

The calibration of your vaporizer and anesthesia system plays a crucial part in your patients health and recovery of surgery. Overtime of the process of evaporation of the anesthetic fluids in the vaporizer result in a varnish-like deposit on the component parts.

Medicanix has the unique ability to check your vaporizers calibration while onsite as well as overhaul the rest of the entire anesthesia system. This service not only saves you money, but also the inconvenience of being without your vaporizer. Should the vaporizer be found to be out-of-spec, we provide a loaner while we take your vaporizer to be overhauled.

We also carry a wide array of anesthesia system parts such as seals, gaskets, pop off valves, and many other parts that may be malfunctioning which are usually shown up during an onsite maintenance.



Microscope Cleaning &
Lubrication Services

Our preventive maintenance service of microscopes includes both its optical and mechanical elements. The microscope's objective lenses and eyepieces are carefully cleaned, the par-focal is adjusted and the prism alignment is checked.

Mechanical Stage parts are cleaned and re-lubricated as needed. Loose fasteners and correct alignment of moving parts are checked and calibrated. The illumination system on the microscope is also checked for efficient operation.

A clean and well lubricated microscope can make all the difference to your veterinary hospital and happier co-workers.



Autoclave/Sterilizer
Maintenance Cleaning

Autoclaves contain lots of moving parts that need periodic preventative maintenance to extend the autoclave's life expectancy. This service is a lot more than just running sterilizer cleaner through your sterilizer our protocol includes:

- Replacement of wear and tear parts such as gaskets and filters
- Chamber & Reservoir Cleaning
- Inspection of electrical/electronic components
- Operational check and calibration

You simply cannot afford the risk of this vital equipment malfunctioning in your busy veterinary hospital.

Autoclave & Sterilizers

Centrifuge

Vaporizers

Surgery Lights

Exam Tables

Doggie Dryers

Oxygen Hoses & Fittings

Cages

Microscopes

Electro-surgery Units

Gurneys

Oxygen Concentrators

Patient Warmers

Otosopes

Scales & Table Scales

Medicanix Inc.

30 Nurney St.

Stamford, CT 06902

Web: www.medicanix.com

Email: sales@medicanix.net

Fax: 203-487-7423

Phone: 800-937-3724 or 203-324-3711

